



Surgical blade

Blade handle for surgical blade

Blade remover

# Surgery

Blade remover S

Blade pliers

Disposable scalpel-Standard

Disposable scalpel-Deep

Disposable scalpel-Mini

Safeshield Scalpel

Safeshield Scalpel vessel

Disposable blood lancet

Disposable blood lancet-SP

Disposable dermatome

Skin graft knife blade

Biopsy punch

Chokkin


Disposable prep razor

MADE IN JAPAN

Surgical blade

Feather surgical blades are available in either stainless steel or non-glare carbon steel that satisfies a wide range of demands of surgeons. The blades are individually wrapped and the packages are gamma sterilized. FEATHER is a leader in technological innovation and the creation of extremely sharp, durable blades.

BLADE CROSS-SECTION COMPARISON

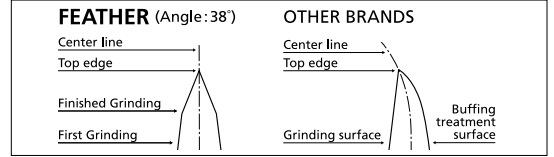

Actual size

	Stainless steel blade	Carbon steel blade
10		
11		
12		
12d		
14		
15		
15s		
15c		
20		
21		
22		
23		
24		
25		

Packaging: Box of 100 pcs.

Actual size

Feather surgical blade No.370,390 and 390C blades are small shaped surgical blades, suitable for microsurgery (fine incision) with microscope. With appropriate shape and thickness, these blades are easy to control and to use for fine incision. Compatible with Feather Surgical Blade Handles (No.3, 3L and 7).

	Stainless steel blade
370	
390	
390C	

Packaging: Box of 10 pcs.

## Blade handle for surgical blade

Various types are available to meet diverse needs. Each is scientifically designed for smooth loading, functional handling ease and perfect balance. The handles are made of stainless steel(except DF-160S-3) to assure long life service and the blade integrates into the handle to form a single instrument. Handles No.3, 7, 3L and DF-160S-3 accommodate blades No.10,11,12,12d,14,15,15c,15s,370,390 and 390C. Handle No.4 and 4L are designed for use with blades No.20,21, 22, 23, 24 and 25.



No.3



No.3 DLC

DLC coating(\*) type is also available.



No.7



No.3L



DF-160S-3 DLC coating(\*)

Round-shaped handle made of titanium alloy. With smooth round shape, great weight balance and ideal length (160mm), this handle provide superior operability for incision in fine/deep area.



No.4



No.4 DLC

DLC coating(\*) type is also available.



No.4L

80% size

\* DLC coating

Feather **DLC coated** blade handles are coated with a diamond like carbon film (DLC) that offers superior performance and durability. The carbon coating prevents the handle from rusting and corroding and the handle design allows for smooth loading of blades.

## Blade remover

Feather Blade remover enables safe and easy removal of used disposable blades from any size or type of scalpel handles. It improves safety and hygiene by eliminating direct contact with blades and facilitates replacement.

Simple - - 1-2-3 operating ease.

Safe & Hygienic - - No direct contact with blades for maximum protection of staff against infection.

Disposable - -When this economical remover becomes full of used blades, just discard it.



● non-sterile

	Blade remover
Dimension(mm)	58(H) x 79.3(W) x 25(D)
Weight(g)	21.7
Blade insertion slot(mm)	17 x 5.2
Material	Polypropylene

● Gamma sterilized  
● Individually packaged  
● Just discard after use.



	Blade remover S
Dimension(mm)	19.3(H)x 79(W)x 27(D)
Weight(g)	11.8
Blade insertion slot(mm)	17 x 5.2
Material	Polypropylene

## Blade pliers

Feather Blade Pliers No.3 enable users to remove a Surgical Blade from the Handle safely. The pliers can be cleaned and resterilized after use.



<CAUTIONS>

- Exclusively for Feather Surgical Blade Handle No.3, 3L, 7 and DF-160S-3. (Cannot be used for No. 4 and 4L.)
- Cannot be used for surgical blade handles of other manufacturers.

	Blade pliers
Weight(g)	40g
Material	Stainless steel
Packaging	1 pc.

**Sterile disposable scalpel**

The Feather disposable scalpel is the last word in precision. Feather, utilizing feedback from surgeons, has developed a design with the sharpest cutting edge available today. Well-balanced construction, remarkable durability and in-the-package sterility further enhance total performance. Best of all, these throw-away scalpels are highly economical.

**Economical & efficient**

Fully automated Feather production and quality control systems are the reasons why these disposable scalpels are the first choice of surgeons everywhere. Feather disposable scalpels provide unsurpassed quality, uniformity and economy.

**Indispensable in emergencies**

Each scalpel is sterile and ready for use as soon as it is removed from the pull-apart aluminum wrapper. You can carry and use them.

**Disposable scalpel - Standard**



70% size

**Disposable scalpel - Deep**

These scalpels feature a plastic handle about twice as long as the standard type for convenient use in deep operations. Available with blades No. 10, 11, 15, 20, 21, 22, 23, 24 and 25.



**Disposable scalpel - Mini**

70% size


This type features a slender plastic handle and is highly suitable for use in surgery in limited areas. For example, the brain, nerves, ears, nose, mouth, plastic surgery, blood vessels, skin and infants. Available with blades No. 11, 14 and 15.



70% size


	Disposable scalpel-Standard	Disposable scalpel-Deep	Disposable scalpel-Mini
10			
11			
12			
12d			
14			
15			
15c			
20			
21			
22			
23			
24			
25			

**How to open the package**

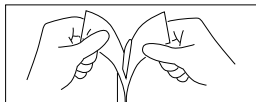

Simply peel the package off by pulling apart the two flaps, as indicated.



**How to remove the snap-off cover**

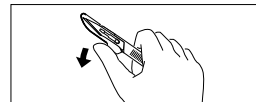



The blade cover is easily removed by applying firm pressure to the snap-tab, as indicated.

## Sterile Safeshield Scalpel


Clear protective shield can be advanced (to cover the blade) and retracted multiple times. 11 sizes are available and color-coded button on the protective shield assists blade identification.

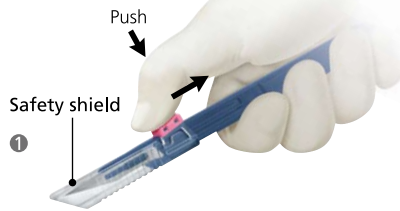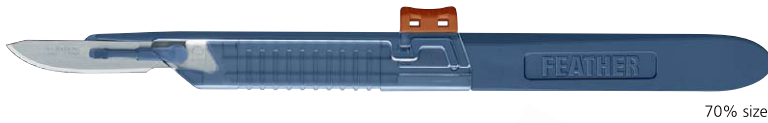

Follow reverse order when disposing of scalpels.

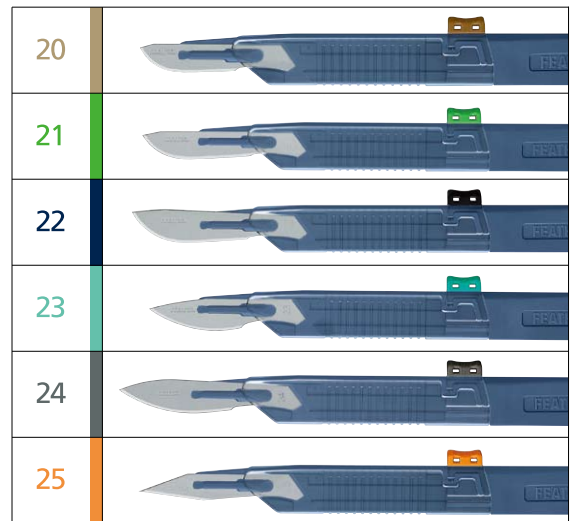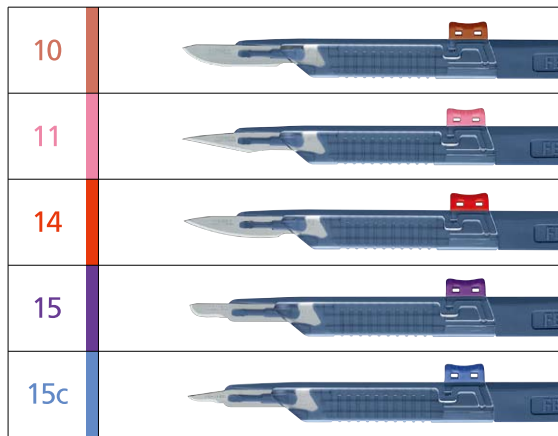



### How to open the package

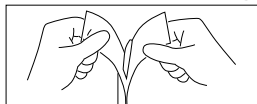

Simply peel the package off by pulling apart the two flaps, as indicated.


## Safeshield Scalpel Vessel

Feather "Safeshield Scalpel Vessel" is designed for custom procedure tray / kit. The safety cover protects the blade in procedure and handling. Supplied in plastic box of 100 pcs, non-sterile.

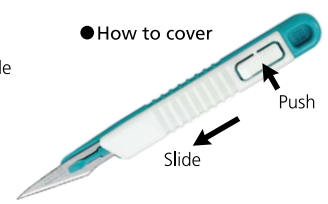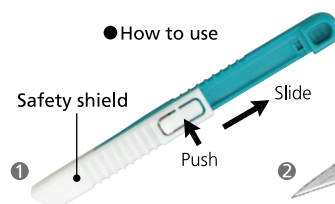

Disposable blood lancet-SP


Feather blood lancet-SP is produced by an original 2 step press treatment that creates an extra sharp cutting edge. The precise, sharp edge allows for smooth operation. Each lancet is packed under strict sanitary conditions. One small box contains 200 lancets (standard packing unit). One export carton contains 40,000 lancets. Each perforated pack contains four lancets.

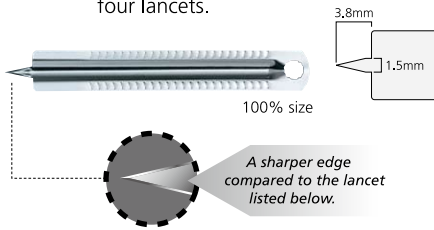

		Disposable blood lancet-SP
Dimension (mm)	40(L) x 5.4(W)	
Material	Stainless steel	
Packaging	One small box contains 200 pcs.	


Disposable blood lancet

Feather Blood Lancet is made of cold-rolled chrome-polished stainless strip steel. Precision engineering assures accuracy to within 1/100mm in thickness, along with superior hardness and point sharpness. Each lancet is packed under strict sanitary conditions. One small box contains 200 lancets (standard packing unit). One export carton contains 40,000 lancets. Each perforated pack contains four lancets.

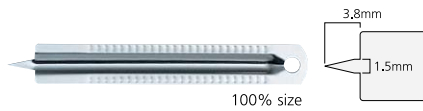

DERMATOME

Disposable dermatome

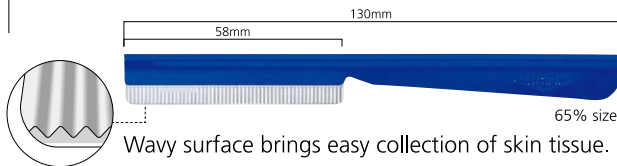

		Disposable dermatome
Dimension	Blade (mm)	58(L) x 16(W)
	Handle (mm)	130(L) x 12(W)
Material	Blade : Stainless steel (specially processed blade)	
Packaging	Sterilized by Gamma radiation at 25kGy. Box of 20 pcs.	

Skin graft knife blade


Feather skin graft knife blade is a disposable blade made of stainless steel. Its sharp edge offers a smooth grafting for the burnt skin. Skin graft knife blade (SGB-157) is compatible with the following knives: ROSENBERG, BRAITHWAITE, COBBETT, and WATSON.



		Skin graft knife blade(SGB-157)
Dimension (mm)	157(L) x 18(W) x 0.229(T)	
Material	Stainless steel	
Packaging	Individually packed in aluminum foil and sterilized by Gamma radiation at 25kGy. Supplied in a box of 20 pcs.	

Biopsy punch

Feather Biopsy Punch offers excellent sharpness by ultra-precision sharpening technology which originates in surgical blade manufacturing.

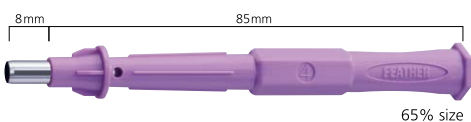

		Biopsy punch
Inner Diameter (mm) / Product No.	1.5 / FBP-15	3.0 / FBP-30
	4.0 / FBP-40	5.0 / FBP-50
Material	Blade	Stainless steel
	Handle	Plastic
Packaging	Sterilized by Gamma radiation at 25 kGy. Box of 20 pcs.	

**Chokkin (with strap)**

Pocket-sized scissors designed for exceptional handling ease, cutting ability and safety. "Chokkin" scissors are small enough to fit in a nurse's pocket. And come with a holder and clip for easy access without damage to the pocket. The perfect tool for the day-to-day hospital routine of nurses. Blade edges are specially processed with fluoroplastic. Residual glue from adhesive tape will not cling to blade after cutting, thus assuring smooth cutting reliability of the scissors.

**Improved Safety**

The blade tip is rounded and protected with a plastic guard. When cutting bandages, only the plastic guard touches the skin. Thus, adding relief and reassurance to the patient during cutting procedures.

**Friendly Design**

"Chokkin" are designed with gentle, curved lines. ("Chokkin" also come in a special "koala" configuration for pediatric nurses that are designed to alleviate anxiety in younger patients.)



\*Stainless steel blade with plastic handles (ABS resin, heatproof temperature : 70°C).

**Disposable prep razor 121 S-B**

- \*Quick and easy to use.
- \*Individually blister packed.
- \*Superbly sharp double-edge blade.



		Disposable prep razor 121S-B	
		Dimension(mm)	65(L)
Material		Blade	Stainless steel
		Handle	Plastic
		Packaging	Supplied in box of 50 pcs.



**FEATHER SAFETY RAZOR CO., LTD.**  
**OVERSEAS TRADE DIVISION**

3-70, OHYODO MINAMI 3-CHOME, KITA-KU, OSAKA 531-0075, JAPAN  
PHONE: +81-6-6458-1638 FAX: +81-6-6458-1611  
URL <http://www.feather.co.jp/> E-mail [overseas@feather.co.jp](mailto:overseas@feather.co.jp)